


2016

Omunis

*We make it **accessible***

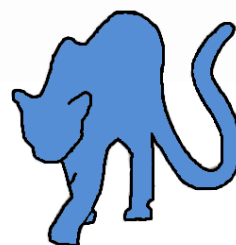

PUMA HCV Kit
PCR. Universal. Molecular. Access.

Cap Alpha – 9 Avenue de l'Europe – 34830 Clapiers, France
Tel : +33 (0)4.34.35.91.18 – www.omunis.com – Email : contact@omunis.com

PUMA HCV Kit



OPEN system

Technical validations

Compared to the *reference technique**

* Commercial PCR, contact us for more details

	OMUNIS	Reference technique
Repeatability	1.8 - 5.8%	1 - 3%
Reproducibility	1.5 - 5.9%	1 - 5%
LOD	80 IU/ml	15 IU/ml
Linearity	$R^2 = 0.9917$ (7log)	$R^2 = 0.9991$
Diversity**	6 genotypes	6 genotypes
Serum volume	140µl	500µl

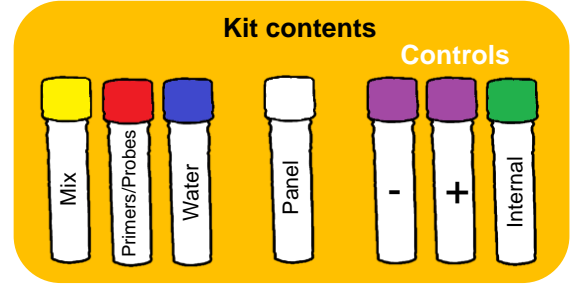
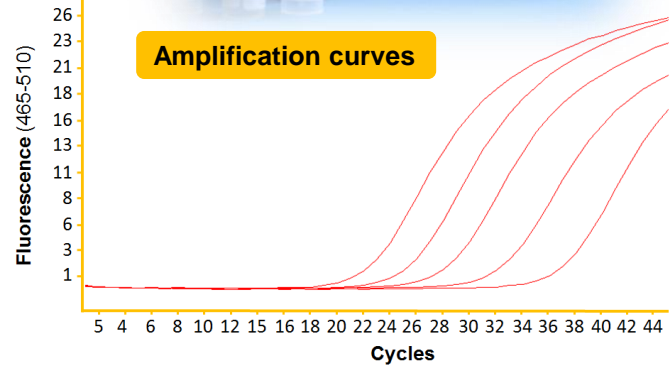


Clinical validations

Compared to the *reference technique**

** (WHO Reference Material 2st HCV RNA Genotype Panel for NAAT)

* Commercial PCR, contact us for more details

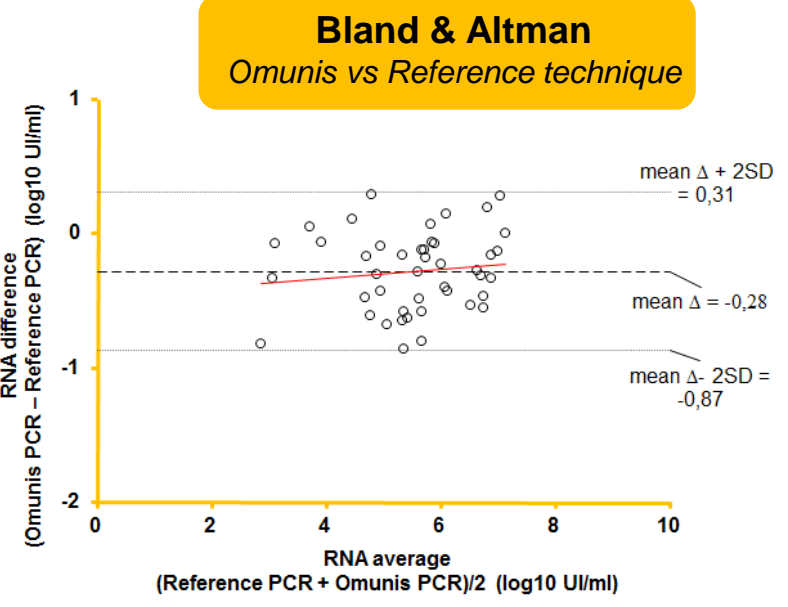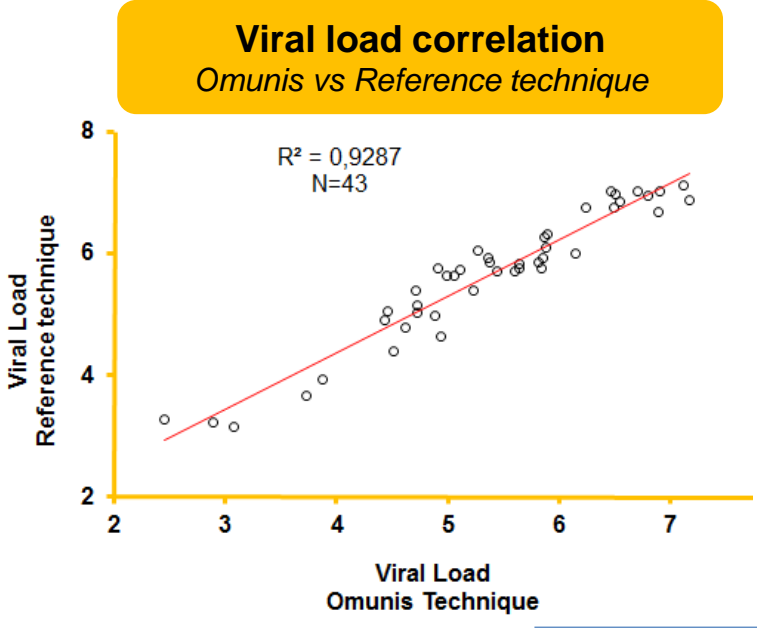

Specificity: 100%
On 64 clinical samples

Sensitivity: 100%
On 64 clinical samples